



Electric Vertical Birthing Chair

Variety of models: SPE / SPE AI / SPE FE AI

Medical technology

SPE R01 VERSION



Features

BENEFITS

The Electric Vertical Birthing Chair is the result of a careful design process, whose development draws from exacting functional requirements. This model offers effective solutions for the widest variety of movements and positions required for delivery. It can perform the standard movements and can be converted into a surgical table when circumstances require it (e.g. for a C-section).

ERGONOMICS

The wide range of movements provides easy access to the operation area, without the need to move the patient. The ergonomic design offers protection and comfort for the patient and easy access to the working area for the surgeon.

MATERIALS

The Electric Birthing Chair is made almost entirely of AISI 304 stainless steel, which provides strength and stability. It has two dual front wheels 125 mm diameter and two rear wheels 125 mm diameter, which enable an agile and smooth movement. It has a central brake which provides the necessary locking function, providing firmness and stability and can accommodate a patient load of up to 200 kg.



Included accessories



Anesthesia screen



Pair of shoulder pads



Pair of grips



Set of 5 clamps



Pair of arm rests



Pair of Goepel leg rests



Pair of leg rests (vertical delivery)



IV pole

VARIETY OF MODELS


AI

304 stainless steel

FE AI


Electric fowler, 304 stainless steel

Movements and positions

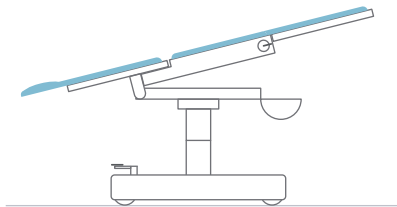

Trendelenburg

14°

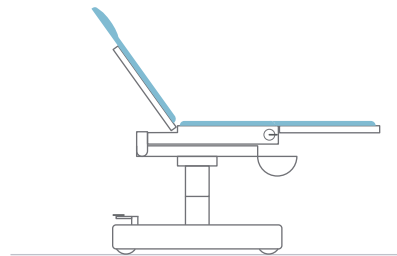


Back rest plate up (SPE and SPE AI models)

47°


Back rest plate up (SPE FE AI model): 55°

Prestaciones y dimensiones

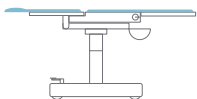

Maximum height

102 cm

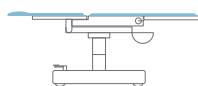

Minimum height

87 cm



Maximum patient load

200 Kg

Chair width	58 cm
Width with side guides	63 cm
Back rest length	77 cm
Seat length	45 cm
Foot panel length	51 cm
Total length	185 cm
Base length	87 cm
Base width	52 cm

DO NOT DISPOSE OF ON PUBLIC ROADS. ND MEDICAL S.R.L. RESERVES THE RIGHT TO MAKE CHANGES WITHOUT PRIOR NOTICE.

Distributed by

Carriego 851, Rosario, Santa Fe, Argentina.
Phone number: +54 341 4361082 / +54 9 341 6206906
info@ndmedical.com.ar

Ndmedical

www.ndmedical.com.ar